




PRESSALIT®

Product Item number	MCT 5 change trolley, manually height adjustable R9502719000
Color	EAN:5708590389875 FDA:INW/D546021 White 
Product description	Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user. MCT 5 change trolley, 27.6" x74.8", two fold-down sides. Manually height adjustable with foot pedal (hydraulic). 19.7" lifting range, from 21,7" - 41.3" Mattress to be ordered separately.
Adjustability Adjustment height Product capacity	Height adjustable 19,7 inch Maximum user weight (MUM): 441 lbs. Maximum load (SWL): 452 lbs. <i>The maximum load is to be understood as the total safe working load of the product, including both the weight of the user, and the weight of the mattress</i> <i>This product has been tested statically with a load of 678 lbs, as the ISO 17966 standard requires testing with a safety factor 1.5.</i>
Materials	Plastics, compact laminate and powder-lacquered steel.

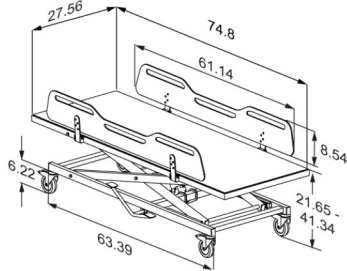

Product	MCT 5 change trolley, manually height adjustable
Item number	R9502719000
Ergonomy	• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user. • One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table. • The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training. • Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.
Features	• The trolley can be moved from room to room for easy transport and multiple users. • Intuitive, easy-to-use drop down rails with one hand operation.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



PRESSALIT®

Product	MCT 5 change trolley, manually height adjustable
Item number	R9502719000
Cleaning	The product is to be cleaned with a damp cloth that has been wrung out in a mild soap mixture, such as universal cleaning solution. After cleaning, it is essential to wipe the product with a dry, soft cloth. This step is crucial to prevent the materials from discoloration. Additionally, it helps prevent any potential degradation of the surface paint, The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättspnit 70% Accel PREvention Wipes Lingettes NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.
Declaration of conformity	This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.

pressalit.com

PRESSALIT®